




Product Specification

Catalog No. 26-45XX-02
5' Digoxigenin-Oligo d(T)

Oligo d(T) 5' Digoxigenin Labeled

For research use only. Not for use in diagnostic procedures for clinical purposes.
Commercial licenses may be obtained directly from Gene Link.

Quantity 25 µg
Shipping Condition Ambient
Storage -20°C

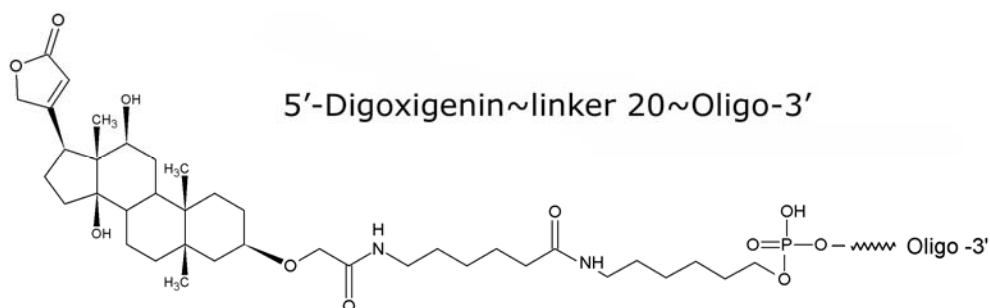

Item	Catalog No.	Product	Sequence (5'-3')	Size	MW*	~nmol/25µg
<input type="checkbox"/>	26-4500-02	5'-Dig-Oligo d(T)12-18	5'-Dig-(T)12-18-3'	12-18	5280	5.2
<input type="checkbox"/>	26-4512-02	5'-Dig-Oligo d(T)12	5'-Dig-(TTTTTTTTTTTT)-3'	12	4367.4	7
<input type="checkbox"/>	26-4513-02	5'-Dig-Oligo d(T)13	5'-Dig-(TTTTTTTTTTTTTT)-3'	13	4671.6	6.5
<input type="checkbox"/>	26-4514-02	5'-Dig-Oligo d(T)14	5'-Dig-(TTTTTTTTTTTTTTT)-3'	14	4975.8	6
<input type="checkbox"/>	26-4515-02	5'-Dig-Oligo d(T)15	5'-Dig-(TTTTTTTTTTTTTTTT)-3'	15	5280	5.5
<input type="checkbox"/>	26-4516-02	5'-Dig-Oligo d(T)16	5'-Dig-(TTTTTTTTTTTTTTTTT)-3'	16	5584.2	5.2
<input type="checkbox"/>	26-4517-02	5'-Dig-Oligo d(T)17	5'-Dig-(TTTTTTTTTTTTTTTTTT)-3'	17	5888.4	5
<input type="checkbox"/>	26-4518-02	5'-Dig-Oligo d(T)18	5'-Dig-(TTTTTTTTTTTTTTTTTTT)-3'	18	6192.6	4.7
<input type="checkbox"/>	26-4519-02	5'-Dig-Oligo d(T)19	5'-Dig-(TTTTTTTTTTTTTTTTTTTT)-3'	19	6496.8	4.4
<input type="checkbox"/>	26-4520-02	5'-Dig-Oligo d(T)20	5'-Dig-(TTTTTTTTTTTTTTTTTTTTT)-3'	20	6801	4.1
<input type="checkbox"/>	26-4521-02	5'-Dig-Oligo d(T)21	5'-Dig-(TTTTTTTTTTTTTTTTTTTTTT)-3'	21	7105.2	3.8

*Molecular weight includes the mw of the dye. An average weight is reported for Oligo d(T)12-18

Description

Oligo d(T)12-18 is the classic primer mix used to prime synthesis of the first strand cDNA by reverse transcriptase using poly A⁺ mRNA as a template. Varying sizes of oligo dT are also used for specific reverse transcription.

Digoxigenin labeled Oligo d(T) primers of varying sizes are individually synthesized with amino C6 linker and conjugated with digoxigenin NHS ester. The digoxigenin conjugated oligos are then gel purified. Digoxigenin labeled Oligo d(T)12-18 are mixed in equimolar ratio. Oligo dT is used primarily to prime synthesis by reverse transcriptase of the first strand cDNA using mRNA as a template for use in microarray hybridization protocol.

The digoxigenin labeled oligo dT is gel purified and supplied as a lyophilized powder, after reconstitution store at -20°C. Oligo purity is greater than 98% as determined by denaturing polyacrylamide gel electrophoresis.

Reconstitution

Recommended reconstitution is at a concentration of ~50 µM in RNase-free DEPC treated water or TE pH 8.

- Spin the tube briefly to bring down the contents of the tube that may have lodged in the cap during shipment. Pellet may be very small and not visible.

- Add 20 X nmol= µl RNase free water or TE directly to the tube. Vortex briefly. See table for exact amount nmol.

Example: for oligo dT 12-18 the nmol supplied is 5.2 (25 µg)

20 X 5.2= 104µl

Add 104 µl RNase free water or TE directly to the tube.

- The above solution is ~50µM. This is equivalent to ~50 pmol/µl and approximately 0.25 µg/µl.

Store at -20°C or below after reconstitution.

Recommended Usage

Use 2 µl of the 50 µM solution for 1 µg poly (A)⁺ RNA as a template in a 20 µl reaction volume. See reaction conditions for more details.

Quality Control Data

This product is certified to prime first strand cDNA reaction using poly (A)⁺ RNA as a template.

Functional Assay Conditions

The conditions given below have been tested to yield first strand cDNA synthesis and is given as an example. Variations and other protocols have been used by other laboratories using this product to yield excellent first strand synthesis. Investigators can substitute their own reaction conditions.

The quality of RNA is very important for the reverse transcription reaction. It is essential to have intact full length RNA as the template material that is free of even trace amounts of RNases and contaminating chemicals. Poor quality RNA template is usually the cause of truncated and incomplete cDNA products.

Add components in the order given below. Reaction volume can be scaled up.

Component	Volume	Comments
poly(A) ⁺ RNA in sterile water Quantity ~1.0 µg	up to 10 µl	Use RNase free reagents and disposables.
RNase-free water	variable	Calculate total volume and add appropriate volume of RNase-free water at this stage.
50 µM oligo(dT)12-18 primer solution (50 pmol/µl = ~0.25 µg/µl)	2 µl	Final concentration is 5 µM (5 pmol/ µl).
Heat mixture to 70°C for 10 min, and quick chill on ice.		
5X first strand buffer [250 mM Tris-HCl (pH 8.3), 375 mM KCl, 15 mM MgCl ₂]	4 µl	
0.1 M DTT	2 µl	
dNTPs (5 mM each dNTP)	2 µl	Final concentration is 0.5 mM of each dNTP.
[α- ³² P]dCTP (1 µCi/µl)	1 µl	Tracer optional. Add only if required.
Reverse transcriptase; 200 units	1- 2 µl	
Total Volume	20 µl	

Incubate at 37°C for 1 hour.

Follow Roche protocol for Digoxigenin based detection system using anti digoxigenin antibody labeled with alkaline phosphatase or other ligands for chemiluminescent or color based protocols.

Related Products

Gene Link stocks various oligo dT primers, oligo dTVN primer, Oligo dT T7 primer, random primers, including an array of fluorescent dye labeled primers for genetic analysis using fluorescent detecting instruments. The C-12 amino labeled primers are ready to be conjugated to the investigators choice of NHS-activated ligand.

Random Primers are a mixture of oligonucleotides representing all possible sequence for that size. Random Primers can be used to prime synthesis in oligo-labeling similar to using hexamers (1,2) and cDNA synthesis. Random prime labeling yields high specific activity labeled DNA probe which can be used for all southern, northern and in situ hybridization studies. Random Primers can be also used similar to using hexamers in cDNA synthesis in combination with oligo dT to yield more 5' end cDNA sequence.

Product Ordering Information

Digoxigenin labeled primers

Catalog No.	Product Description	Quantity
26-4500-02	5'-Dig-Oligo d(T)12-18	25ug
26-4518-02	5'-Dig-Oligo d(T)18	25ug
26-4519-02	5'-Dig-Oligo d(T)19	25ug
26-4520-02	5'-Dig-Oligo d(T)20	25ug
26-4521-02	5'-Dig-Oligo d(T)21	25ug

Related Product Ordering Information

Oligo dT unlabeled primers

Catalog No.	Product Description	Quantity
26-4000-04	Oligo d(T) 12	100 µg
26-4000-01	Oligo d(T)16	100 µg
26-4000-02	Oligo d(T)18	100 µg
26-4000-05	Oligo d(T)12-18	100 µg
26-4002-10	Oligo d(T)23	50 µg
26-4002-11	Oligo d(T)23 VN	50 µg
26-4002-16	Oligo d(T)36	50 µg
26-3000-23	T7 Oligo d(T)23	50 µg
26-3000-24	T7 Oligo d(T) 23 VN	50 µg

Related Product Ordering Information

Oligo dT Cy5 labeled primers

Catalog No.	Product Description	Quantity
26-4400-02	5'-Cy5 Oligo d(T)12-18	25ug
26-4412-02	5'-Cy5 Oligo d(T)12	25ug
26-4413-02	5'-Cy5 Oligo d(T)13	25ug
26-4414-02	5'-Cy5 Oligo d(T)14	25ug
26-4415-02	5'-Cy5 Oligo d(T)15	25ug
26-4416-02	5'-Cy5 Oligo d(T)16	25ug
26-4417-02	5'-Cy5 Oligo d(T)17	25ug
26-4418-02	5'-Cy5 Oligo d(T)18	25ug
26-4419-02	5'-Cy5 Oligo d(T)19	25ug
26-4420-02	5'-Cy5 Oligo d(T)20	25ug
26-4421-02	5'-Cy5 Oligo d(T)21	25ug

Visit www.genelink.com for a complete list of fluorescent dye labeled oligo dT primers

Related Product Ordering Information

Random Primers

Catalog No.	Product Description	Quantity
26-4000-03	Random Hexamers	100ug
26-4000-06	Random Nonamers	100ug
26-4000-07	Random Heptamer Phosphorylated pd(N)7	50ug
26-4000-08	Random Octamer Phosphorylated pd(N)8	50ug
26-4000-09	Random Nonamer Phosphorylated pd(N)9	50ug
26-4000-10	Random Hexamer Phosphorylated pd(N)6	50ug
26-4000-11	Random Heptamer	100ug
26-4000-12	Random Octamer	100ug
26-4000-13	Random 12mers	100ug
26-4000-16	Random 15mer	100ug
26-4000-14	Random 24mers	100ug
26-4000-15	Random 36mers	100ug
26-4000-17	Random 60mer	100ug
26-4001-13	Random Hexamer 72%GC	100ug
26-4001-16	Random Nonamers 72%GC	100ug
26-4001-17	Random 36mer 72%GC	100ug
26-4001-18	Random 60mer 72%GC	100ug
26-4000-81	5'-Dig Random Hexamer	25ug
26-4000-82	5'-Dig Random Heptamer	25ug
26-4000-83	5'-Dig Random Octamer	25ug
26-4000-84	5'-Dig Random Nonamer	25ug
26-4000-91	5'-Amino C12 Random Hexamer	25ug
26-4000-92	5'-Amino C12 Random Heptamer	25ug
26-4000-93	5'-Amino C12 Random Octamer	25ug
26-4000-94	5'-Amino C12 Random Nonamer	25ug
26-4001-01	5'-Biotin Random Hexamer	25ug
26-4001-02	5'-Biotin Random Heptamer	25ug
26-4001-03	5'-Biotin Random Octamer	25ug
26-4001-04	5'-Biotin Random Nonamer	25ug

Random Primers Product Ordering Information

Catalog No.	Product Description	Size
26-4000-21	5'-Cy3 Random Hexamer	25ug
26-4000-22	5'-Cy3 Random Heptamer	25ug
26-4000-23	5'-Cy3 Random Octamer	25ug
26-4000-24	5'-Cy3 Random Nonamer	25ug
26-4000-26	5'-Cy3 Random 36mer	25ug
26-4000-25	5'-Cy3 Random 60mer	25ug
26-4001-23	5'-Cy3 Random Hexamers 72%GC	25ug
26-4001-26	5'-Cy3 Random Nonamers 72%GC	25ug
26-4001-27	5'-Cy3 Random 36mers 72%GC	25ug
26-4001-28	5'-Cy3 Random 60mers 72%GC	25ug
26-4000-31	5'-Cy5 Random Hexamer	25ug
26-4000-32	5'-Cy5 Random Heptamer	25ug
26-4000-33	5'-Cy5 Random Octamer	25ug
26-4000-34	5'-Cy5 Random Nonamer	25ug
26-4000-36	5'-Cy5 Random 36mer	25ug
26-4000-35	5'-Cy5 Random 60mer	25ug
26-4001-33	5'-Cy5 Random Hexamers 72%GC	25ug
26-4001-36	5'-Cy5 Random Nonamers 72%GC	25ug
26-4001-37	5'-Cy5 Random 36mers 72%GC	25ug
26-4001-38	5'-Cy5 Random 60mers 72%GC	25ug
26-4000-41	5'-HEX Random Hexamer	25ug
26-4000-42	5'-HEX Random Heptamer	25ug
26-4000-43	5'-HEX Random Octamer	25ug
26-4000-44	5'-HEX Random Nonamer	25ug
26-4000-51	5'-FAM Random Hexamer	25ug
26-4000-52	5'-FAM Random Heptamer	25ug
26-4000-53	5'-FAM Random Octamer	25ug
26-4000-54	5'-FAM Random Nonamer	25ug
26-4000-61	5'-TET Random Hexamer	25ug
26-4000-62	5'-TET Random Heptamer	25ug
26-4000-63	5'-TET Random Octamer	25ug
26-4000-64	5'-TET Random Nonamer	25ug
26-4000-71	5'-FI Random Hexamer	25ug
26-4000-72	5'-FI Random Heptamer	25ug
26-4000-73	5'-FI Random Octamer	25ug
26-4000-74	5'-FI Random Nonamer	25ug