



## Product Specifications

Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

## Oligo Modifications


For research use only. Not for use in diagnostic procedures for clinical purposes.

### Fam-NHS (6-Fam NHS)

Category	Fluorescent Dyes
Modification Code	Fam-N
Reference Catalog Number	26-6730
5 Prime	Y
3 Prime	Y
Internal	Y
Molecular Weight(mw)	358

Fam NHS  
[26-6730-XX]

Amino Linker

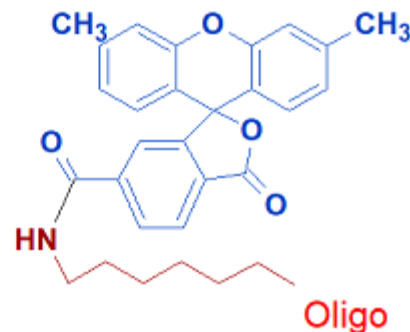

[Click here for a list of fluorophores.](#)

This modification is a post synthesis conjugation to a primary amino group thus an additional modification with an amino group is required. A C3, C6 or C12 amino group can be placed at the 5' or for the 3' end a C3 or C7 amino and for internal positions an amino modified base is used, e.g Amino dT C6.

Yield of Post Synthesis NHS, Maleimide & Click Ligand Conjugation\* Oligo Scale of Synthesis Yield, nmols 50 nmol 2 nmol 200 nmol 5 nmol 1 umol 16 nmol 2 umol 30 nmol 5 umol 75 nmol 10 umol 150 nmol 15 umol 225 nmol \* The yield will be lower for oligos longer than 50mer. [Click here for yield table of long oligos.](#) \* [Click here for RNA Oligos scale of synthesis and yield.](#) **NHS Ligand conjugation** requires a primary amino group. Gene Link offers a wide selection of amino modifications for 5', 3' and internal sites. [Click here for a list of conjugation chemistry modifications.](#) **Maleimide Ligand conjugation** requires a thiol group. Gene Link offers a wide selection of thiol modifications for 5', 3' and internal sites.

Click here for a list of conjugation chemistry modifications. **Click Chemistry Ligand conjugation** requires a corresponding Click modification; examples Alkyne:Azide, Azide:DBCO, BCN:Azide, BCN: TCO:Tetrazine. Gene Link offers a wide selection of click modifications for 5', 3' and internal sites. Click here for a list of click chemistry modifications.

6-carboxyfluorescein (6-FAM) is the most commonly used fluorescent dye for labeling oligonucleotides. 6-FAM is reactive, water-soluble, and has an absorbance maximum of 492 nm and an emission maximum of 517 nm. 6-FAM plays a particularly important role in real-time PCR applications, being used as a reporter moiety in TaqMan probes (1), Scorpion primers (2) and Molecular Beacons (3). For such probes, 6-FAM is most commonly paired with the dark quencher BHQ-1, as the two have excellent spectral overlap. 6-FAM-labeled primers have also been used for bacterial SNP genotyping by allele-specific real-time PCR (4).

6-FAM can be used to label DNA oligos for use as hybridization probes in a variety of in vivo and in vitro research or diagnostic applications, as well as for structure-function studies of DNA, RNA, and protein-oligonucleotide complexes. Oligos labeled with 6-FAM at the 5'-end can be used as PCR and DNA sequencing primers to generate fluorescently-labeled PCR, sequencing or genetic analysis (AFLP or microsatellite) products.

## References

1. Livak, K.J., Flood, S.J.A., Marmaro, J., Giusti, W., Deetz, K. Oligonucleotides with fluorescent dyes at opposite ends provide a quenched probe system useful for detecting PCR product and nucleic acid hybridization. *PCR Methods Appl.* (1995), **4**: 1-6.
2. Thelwell, N., Millington, S., Solinas, A., Booth, J., Brown, T. Mode of action and application of Scorpion primers to mutation detection. *Nucleic Acids Res.* (2000), **28**: 3752-3761.
3. Tyagi, S., Kramer, F.R. Molecular beacons: probes that fluoresce upon hybridization. *Nat. Biotechnol.* (1996), **14**: 303-308.
4. Huygens, F., Inman-Bamber, J., Nimmo-G.R., Munckhof, W., Schooneveldt, J., Harrison, B., McMahon, J.A., Giffard, P.M. Staphylococcus aureus Genotyping Using Novel Real-Time PCR Formats. *J. Clin. Microbiol.* (2006), **44**: 3712-3718.