



Product Specifications


Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

Oligo Modifications

For research use only. Not for use in diagnostic procedures for clinical purposes.

HEX-SIMA dT

Category	Fluorescent Dyes
Modification Code	HEX-SIMA-dT
Reference Catalog Number	26-6499
5 Prime	Y
3 Prime	Y
Internal	Y
Molecular Weight(mw)	1037.79

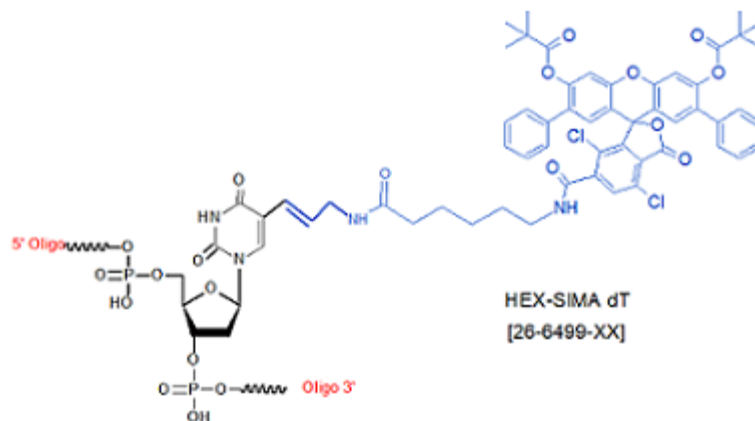

Click here for a list of fluorophores.

HEX-SIMA is Dichloro-diphenyl-fluorescein and exhibits virtually identical absorbance and emission spectra to HEX. HEX-SIMA dT can be placed at any site of an oligo. This is specially suited for internal positions. Hexachloro-fluorescein (HEX) is hexachlorinated version of the fluorescent dye fluorescein, and is used for labeling oligonucleotides at either the 5' or 3' end. HEX has an absorbance maximum of 535 nm and an emission maximum of 556 nm. HEX can be used in real-time PCR applications as a reporter moiety in TaqMan probes (1), Scorpion primers (2) and Molecular Beacons (3). For such probes, HEX is most commonly paired with the dark quencher BHQ-1, as the two have good spectral overlap.

HEX-HEX also can be used to label DNA oligos for use as hybridization probes in a variety of in vivo and in vitro research or diagnostic applications, as well as for structure-function studies of DNA, RNA, and protein-oligonucleotide complexes. Oligos labeled with HEX at the 5' end can be used as PCR and DNA sequencing primers to generate fluorescently-labeled PCR, sequencing or genetic analysis (AFLP or microsatellite) products.

Applied Biosystems Proprietary Dyes & Possible Substitutions

Dye

Color

Absorbance max (nm)

Emission max (nm) VIC Pink Red 538 554 Cal Orange 560 Pink Red 537 558 HEX Pink Red 535 556 NED Red Orange 546 575 Cy3 Red Orange 550 570 PET Red Orange 558 595 Cy3.

5 Red 588 604 ROX Red 575 602 Texas Red Red 583 603

Click here for a list of fluorophores.

References

1. Livak, K.J., Flood, S.J.A., Marmaro, J., Giusti, W., Deetz, K. Oligonucleotides with fluorescent dyes at opposite ends provide a quenched probe system useful for detecting PCR product and nucleic acid hybridization. *PCR Methods Appl.* (1995), **4**: 1-6.
2. Thelwell, N., Millington, S., Solinas, A., Booth, J., Brown, T. Mode of action and application of Scorpion primers to mutation detection. *Nucleic Acids Res.* (2000), **28**: 3752-3761.
3. Tyagi, S., Kramer, F.R. Molecular beacons: probes that fluoresce upon hybridization. *Nat. Biotechnol.* (1996), **14**: 303-308.