



Product Specifications


Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

Oligo Modifications

For research use only. Not for use in diagnostic procedures for clinical purposes.

HEX-N

Category	Fluorescent Dyes
Modification Code	Hex-N
Reference Catalog Number	26-6590
5 Prime	Y
3 Prime	Y
Internal	Y
Molecular Weight(mw)	744.13

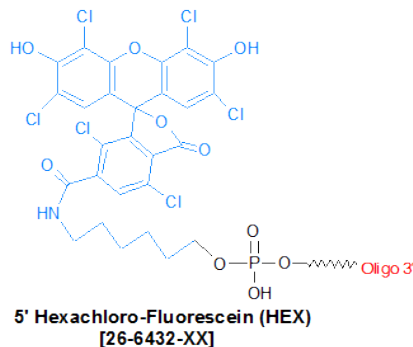

This modification is a post synthesis conjugation to a primary amino group thus an additional modification with an amino group is required. A C3, C6 or C12 amino group can be placed at the 5' or for the 3' end a C3 or C7 amino and for internal positions an amino modified base is used, e.g Amino dT C6.

Yield of Post Synthesis NHS, Maleimide & Click Ligand Conjugation* Oligo Scale of Synthesis Yield, nmols 50 nmol 2 nmol 200 nmol 5 nmol 1 umol 16 nmol 2 umol 30 nmol 5 umol 75 nmol 10 umol 150 nmol 15 umol 225 nmol * The yield will be lower for oligos longer than 50mer. Click here for yield table of long oligos. * Click here for RNA Oligos scale of synthesis and yield. **NHS Ligand conjugation** requires a primary amino group. Gene Link offers a wide selection of amino modifications for 5', 3' and internal sites. Click here for a list of conjugation chemistry modifications. **Maleimide Ligand conjugation** requires a thiol group. Gene Link offers a wide selection of thiol modifications for 5', 3' and internal sites.

Click here for a list of conjugation chemistry modifications. **Click Chemistry Ligand conjugation** requires a corresponding Click modification; examples Alkyne:Azide, Azide:DBCO, BCN:Azide, BCN: TCO:Tetrazine. Gene Link offers a wide selection of click modifications for 5', 3' and internal sites. Click here for a list of click chemistry modifications.

Hexachloro-fluorescein (HEX) is hexachlorinated version of the fluorescent dye fluorescein, and is used for labeling oligonucleotides at either the 5' or 3' end. HEX has an absorbance maximum of 535 nm and an emission maximum of 556 nm. HEX can be used in real-time PCR applications as a reporter moiety in TaqMan probes (1), Scorpion primers (2) and Molecular Beacons (3). For such probes, HEX is most commonly paired with the dark quencher BHQ-1, as the two have good spectral overlap.

HEX also can be used to label DNA oligos for use as hybridization probes in a variety of in vivo and in vitro research or diagnostic applications, as well as for structure-function studies of DNA, RNA, and protein-oligonucleotide complexes. Oligos labeled with HEX at the 5' end can be used as PCR and DNA sequencing primers to generate fluorescently-labeled PCR, sequencing or genetic analysis (AFLP or microsatellite) products. NOTE: If HEX is on the 3' end of the oligo, it cannot be used as a primer in PCR-based applications. **References**

1. Livak, K.J., Flood, S.J.A., Marmaro, J., Giusti, W., Deetz, K. Oligonucleotides with fluorescent dyes at opposite ends provide a quenched probe system useful for detecting PCR product and nucleic acid hybridization. *PCR Methods Appl.* (1995), **4**: 1-6.
2. Thelwell, N., Millington, S., Solinas, A., Booth, J., Brown, T. Mode of action and application of Scorpion primers to mutation detection. *Nucleic Acids Res.* (2000), **28**: 3752-3761.
3. Tyagi, S., Kramer, F.R. Molecular beacons: probes that fluoresce upon hybridization. *Nat. Biotechnol.* (1996), **14**: 303-308.