



## Product Specifications


Custom Oligo Synthesis, antisense oligos, RNA oligos, chimeric oligos, Fluorescent dyes, Affinity Ligands, Spacers & Linkers, Duplex Stabilizers, Minor bases, labeled oligos, Molecular Beacons, siRNA, phosphonates Locked Nucleic Acids (LNA); 2'-5' linked Oligos

## Oligo Modifications

For research use only. Not for use in diagnostic procedures for clinical purposes.

### Azide PEG3 Maleimide

Category	Click Chemistry
Modification Code	N3-PEG3-Mal
Reference Catalog Number	26-6761
5 Prime	Y
3 Prime	Y
Internal	Y
Molecular Weight(mw)	369.37

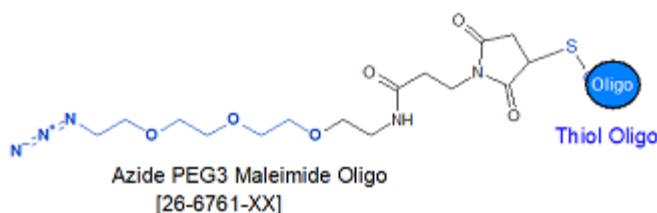

#### [Click here for a complete list of Click Chemistry Oligo Modifications](#)

This modification is a post synthesis maleimide conjugation to a reduced thiol amino group thus an additional modification with thiol group is required. A C3 or C6 thiol group can be placed at the 5' or for internal positions Thiol C6 dT modified base is used. Yield of Post Synthesis NHS, Maleimide & Click Ligand Conjugation\* Oligo Scale of Synthesis Yield, nmols 50 nmol 2 nmol 200 nmol 5 nmol 1 umol 16 nmol 2 umol 30 nmol 5 umol 75 nmol 10 umol 150 nmol 15 umol 225 nmol \* The yield will be lower for oligos longer than 50mer. [Click here for yield table of long oligos.](#) \* [Click here for RNA Oligos scale of synthesis and yield.](#) **NHS Ligand conjugation** requires a primary amino group. Gene Link offers a wide selection of amino modifications for 5', 3' and internal sites.

[Click here for a list of conjugation chemistry modifications.](#) **Maleimide Ligand conjugation** requires a thiol group. Gene Link offers a wide selection of thiol modifications for 5', 3' and internal sites.

Click here for a list of conjugation chemistry modifications. **Click Chemistry Ligand conjugation** requires a corresponding Click modification; examples Alkyne:Azide, Azide:DBCO, BCN:Azide, BCN:Tetrazine and TCO:Tetrazine. Gene Link offers a wide selection of click modifications for 5', 3' and internal sites. Click here for a list of click chemistry modifications.

Azide PEG3 Maleimide can be used to introduce an active azide group to a thiol-modified oligonucleotide. The Azide C2 PEG3 Maleimide is then manually attached to the oligo through the thiol group in a separate reaction post-synthesis. The presence of the azide allows the user to use "Click Chemistry" (a [3+2] cycloaddition reaction between alkynes and azides, using copper (I) iodide as a catalyst) to conjugate the azide-modified oligo to a terminal alkyne-modified oligo with extremely high regioselectivity and efficiency (1,2). Preparation of the alkyne-modified oligo can be achieved using the 5'Hexynyl modifier (see its respective tech sheet for details). Click chemistry can be used to form short, cyclic oligos that can be used as research tools in various biophysical and biological studies (3). In particular, they have considerable potential for in vivo work, as cyclic oligos are known to be very stable in serum for up to several days.

#### References

1. Huisgen, R. *Angew. Chem. Int. Ed.* (1963), **2**: 565-568.
2. Rostovtsev, V.V., Green, L.G., Fokin, V.V., Sharpless, K.B. A Stepwise Huisgen Cycloaddition Process: Copper(I)-Catalyzed Regioselective Ligation of Azides and Terminal Alkynes. *Angew. Chem. Int. Ed.* (2002), **41**: 2596-2599.
3. Kumar, R., El-Sagheer, A., Tumpene, J., Lincoln, P., Wilhelmsson, L.M., Brown, T. Template-Directed Oligonucleotide Strand Ligation, Covalent Intramolecular DNA Circularization and Catenation Using Click Chemistry. *J. Am. Chem. Soc.* (2007), **129**: 6859-6864.